

Launching Kaluza Analysis Software



The shortcut for Kaluza software was created on your desktop during the installation process (as described in the software CD package). To launch the software, double-click the Kaluza icon.

Using the License Key

A license key enables the use of Kaluza after the trial period has ended. License keys are provided by HASP®.

Setting-Up a Computer Using a Single License Key

To set up a single license key:

- 1 Install Kaluza on your computer using the instructions included in the software CD package.
- 2 Plug the USB key into host computer USB port. This allows full access to Kaluza.

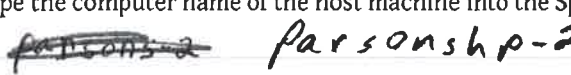
NOTE For additional information regarding the HASP key, refer to the instructions on the website at <http://localhost:1947>.


Setting-Up Computers Using a Network License Key

Prior to setting up user computers on a network license, you must set up the host computer. Follow the instructions in [CHAPTER 1, Setting-Up a Computer Using a Single License Key](#), to complete host computer setup.

To connect network computers to the host computer:

- 1 Install Kaluza on all computers that need to run Kaluza.
- 2 Open a web browser program.
- 3 Enter the following address into the address bar: **<http://localhost:1947>**
You are now connected to the HASP License Manager Admin Control Center.
- 4 From the Administration Options section, select **Configuration**.

- 5 Select the **Access to Remote License Managers** tab.
- 6 Select the **Allow Access to Remote Licenses** check box.
- 7 Type the computer name of the host machine into the **Specify Search Parameters** field.

- 8 Select **Submit**, which connects the computer to the network license key and grants full access to Kaluza.

NOTE To verify that a computer is connected to the network license key, select the  icon after launching Kaluza; this initiates the About screen. In the **License Type** section of the screen, a **Network** license type is indicated when the network license key is recognized by the computer.

NOTE For additional information regarding the HASP key, refer to the instructions on the website at <http://localhost:1947>.

License Key Troubleshooting

If you currently have a HASP license key but are unable to access Kaluza due to a license expiration error similar to the one shown in [Figure 1.1](#), your computer's virus scanner may be preventing access to the HASP license service. To enable access, contact your local Technical Support personnel to request to permission for **hasplms** (HASP License Manager) service on your computer.

Figure 1.1 Kaluza License Error Message

